

SCALPELS AND KNIVES



MICRO SCISSORS



SCISSORS



FORCEPS – ANATOMICS, SURGICALS, MICRO



ANEURYSM VESSEL CLIPS



ARTERY FORCEPS, DRESSING AND SPONGE FORCEPS

RETRACTORS, SELF-RETAINING RETRACTORS, ABDOMINAL RETRACTORS



PROBES, ASPIRATION





SUTURE



DRESSING, PLASTER



BONE SURGERY



BONE SURGERY



THORAX, CARDIOVASCULAR SURGERY



NEUROSURGERY, LAMINECTOMY



STOMACH. INTESTINES, UROLOGY



GYNAECOLOGY, OBSTETRICS



OPHTHALMOLOGY




OTOLOGY



RHINOLOGY



ORAL MAXILLO-FACIAL SURGERY

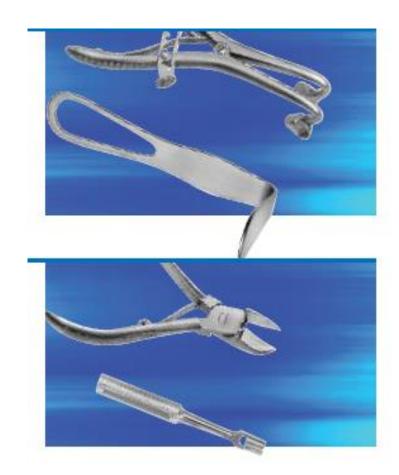


DERMATOLOGY




STERILISATION, HYGIENE, STORAGE

FORCEPS

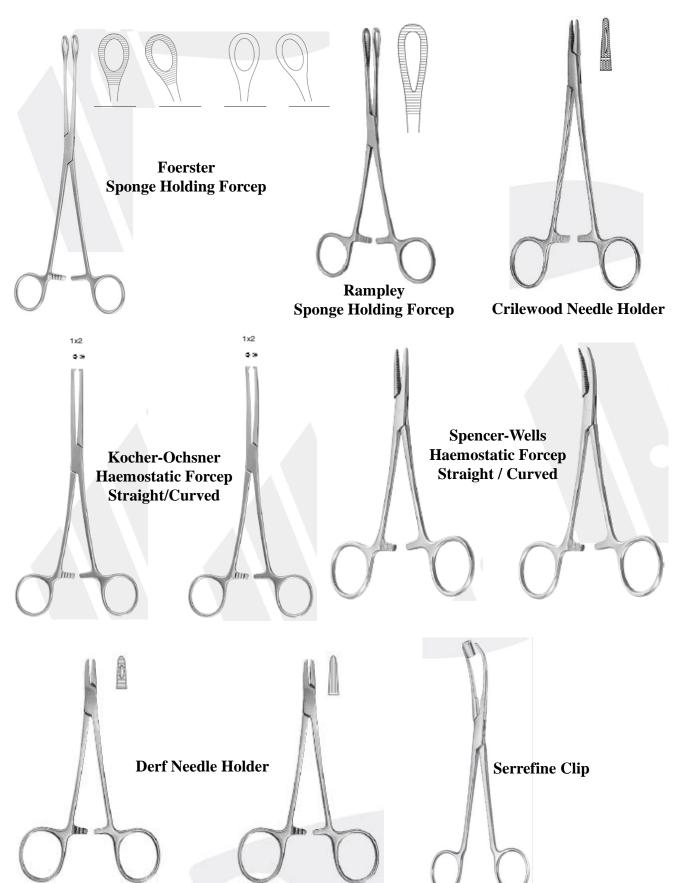


DISSECTING FORCEPS

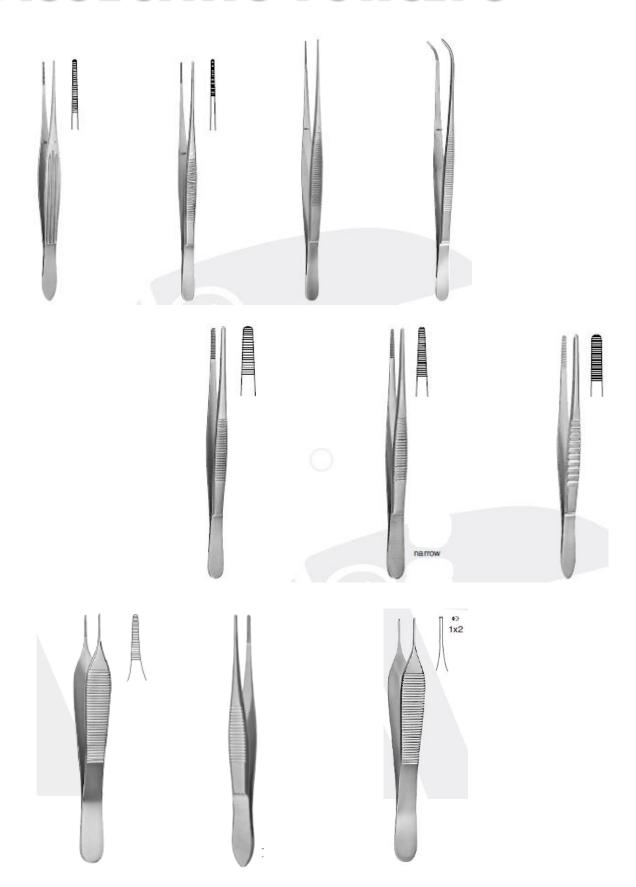

TISSUE FORCEPS WITH TOOTHING




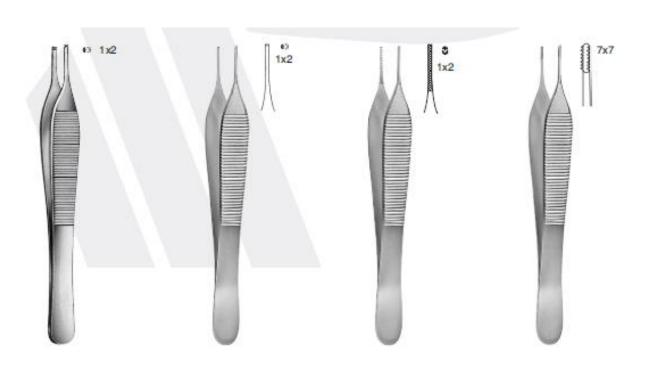


MICRO FORCEPS




SURGICAL, DISSECTING SCISSORS

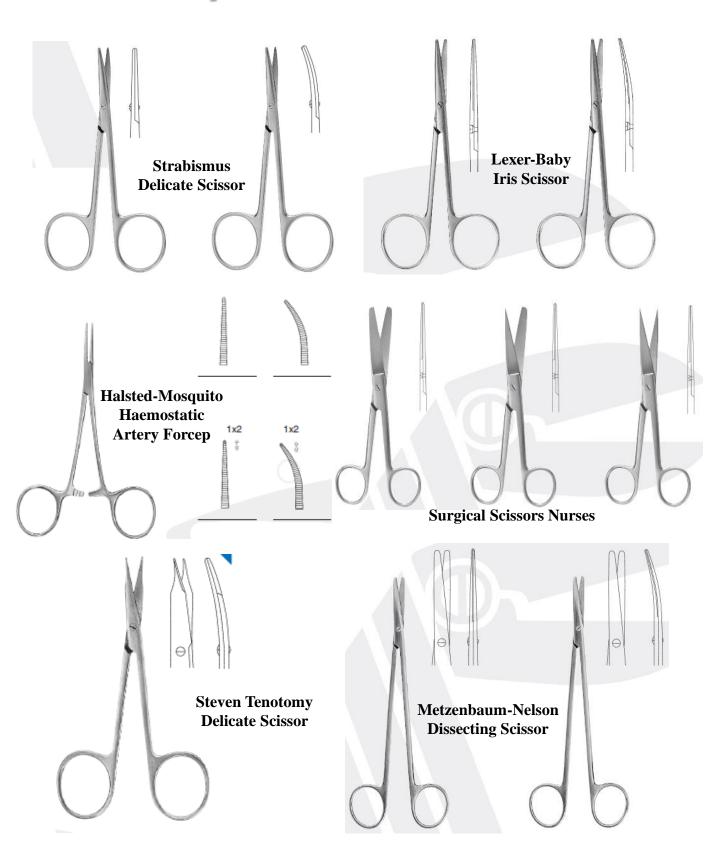


SURGICAL, DISSECTING SCISSORS

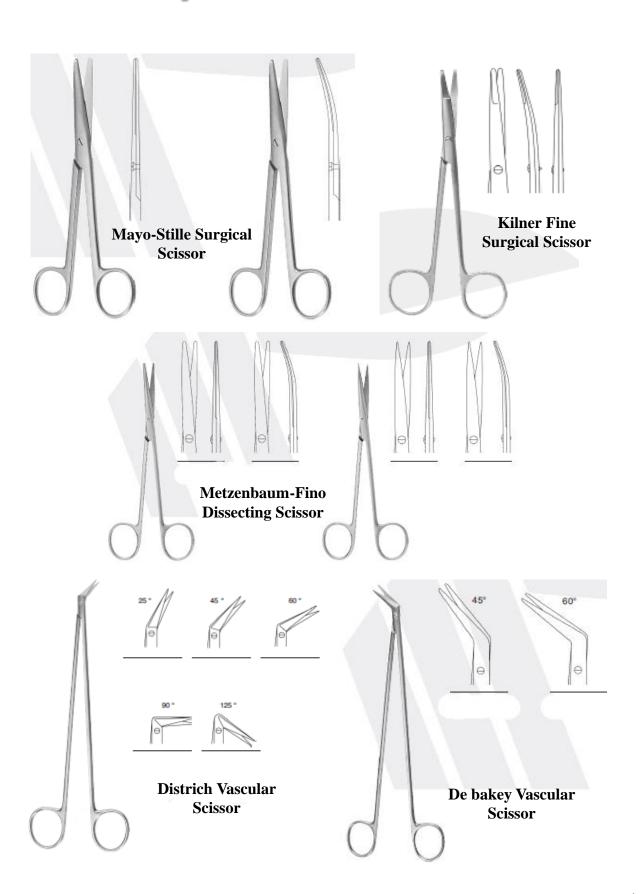


MICRO SCISSORS, NEEDLE HOLDERS

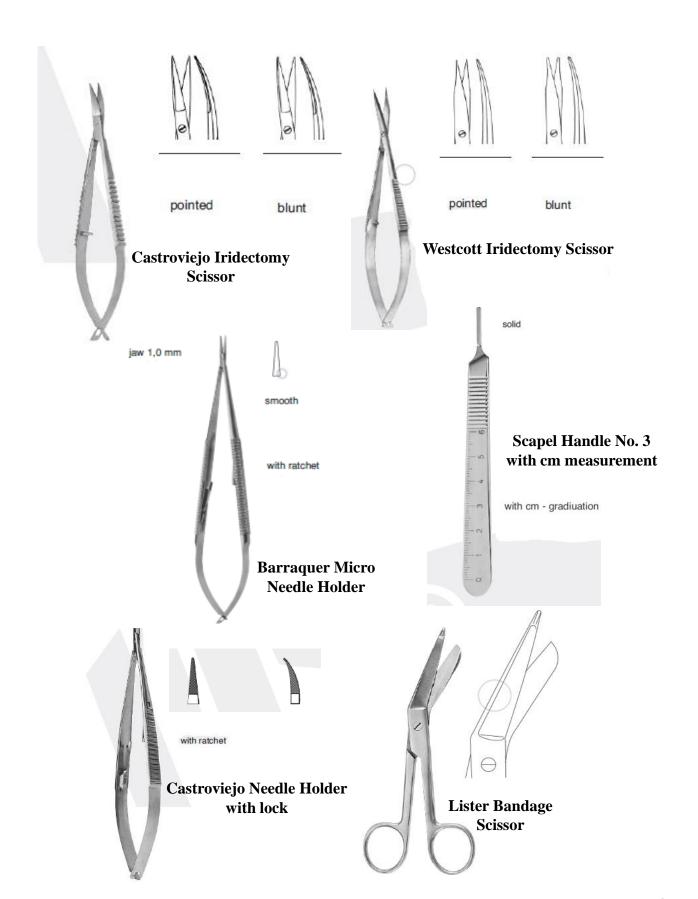

TUNGSTEN CARBIDE INSTRUMENTS

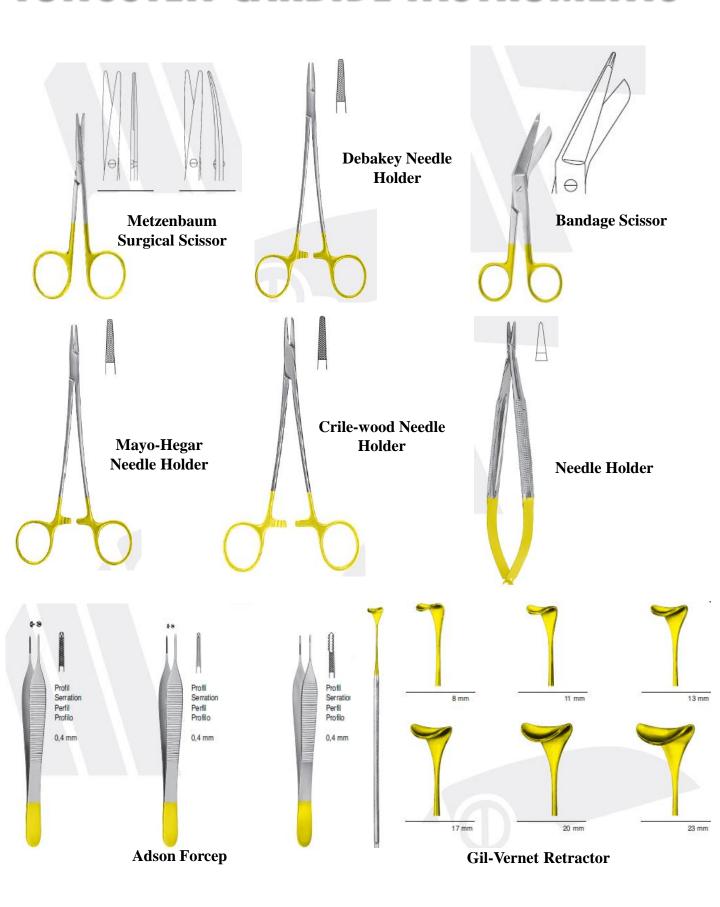


SUPER-CUT SCISSORS

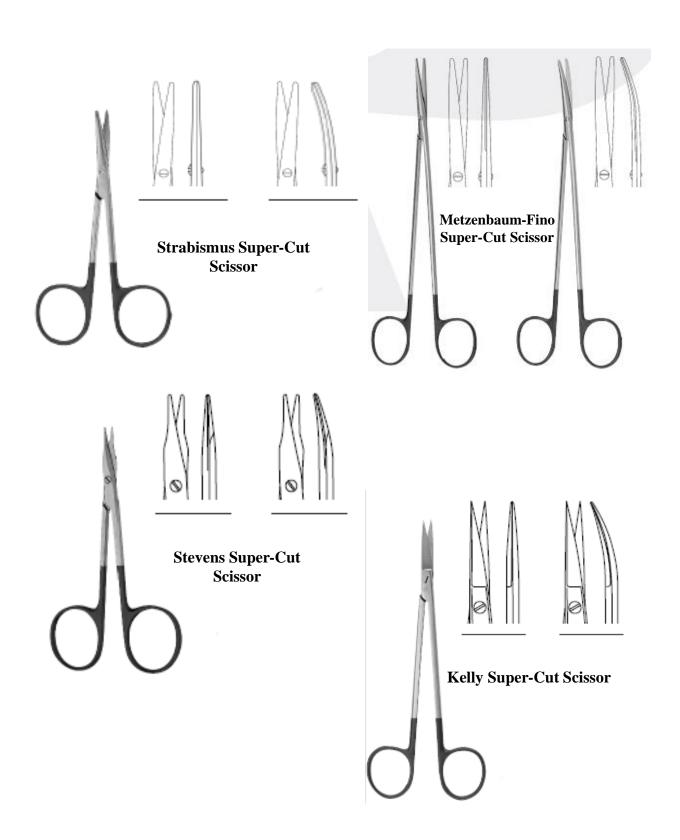


CLAMPS




RETRACTORS

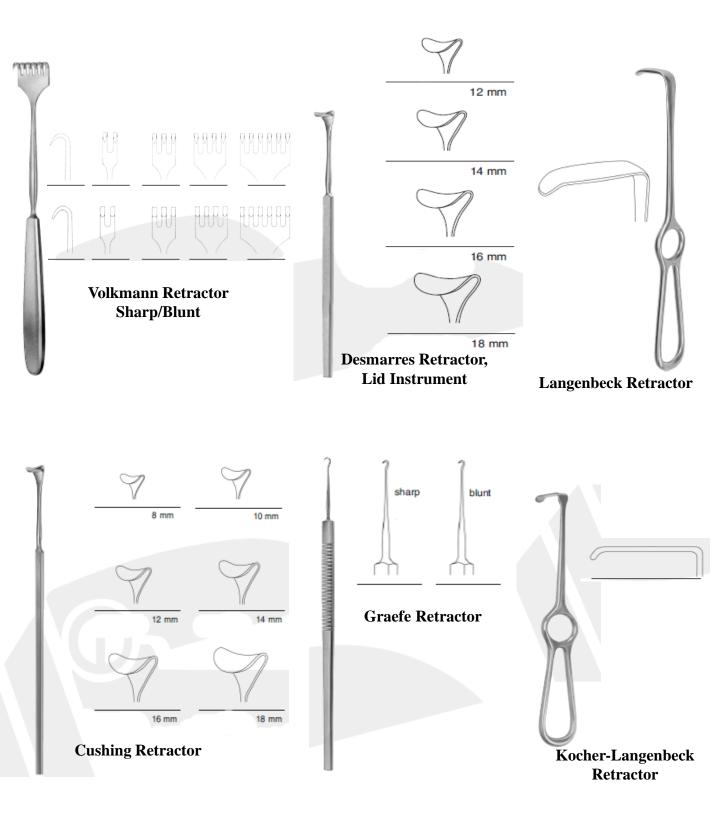



RETRACTORS

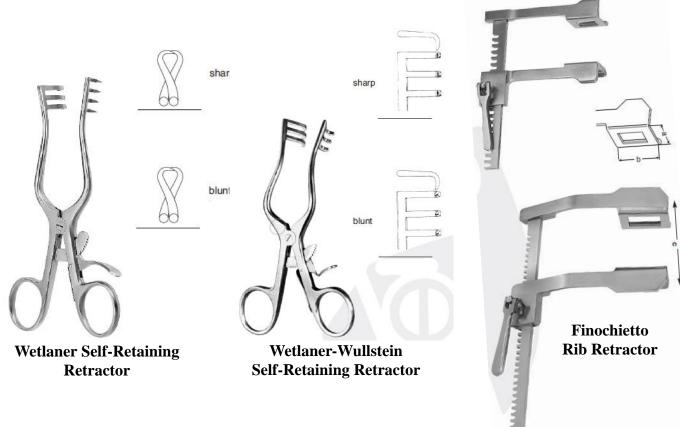

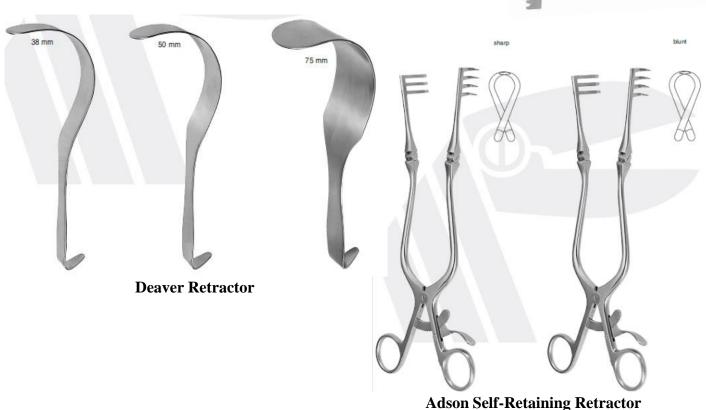


BONE RONGEURS, CUTTING FORCEPS

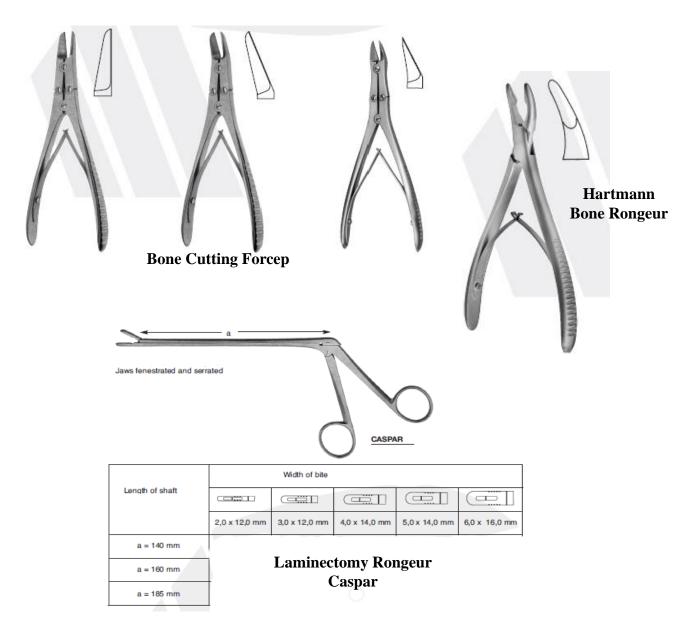


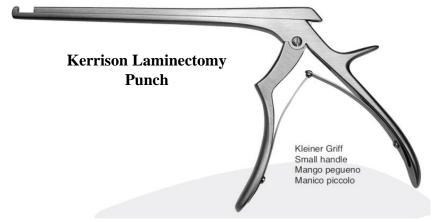



MICRO EAR FORCEPS

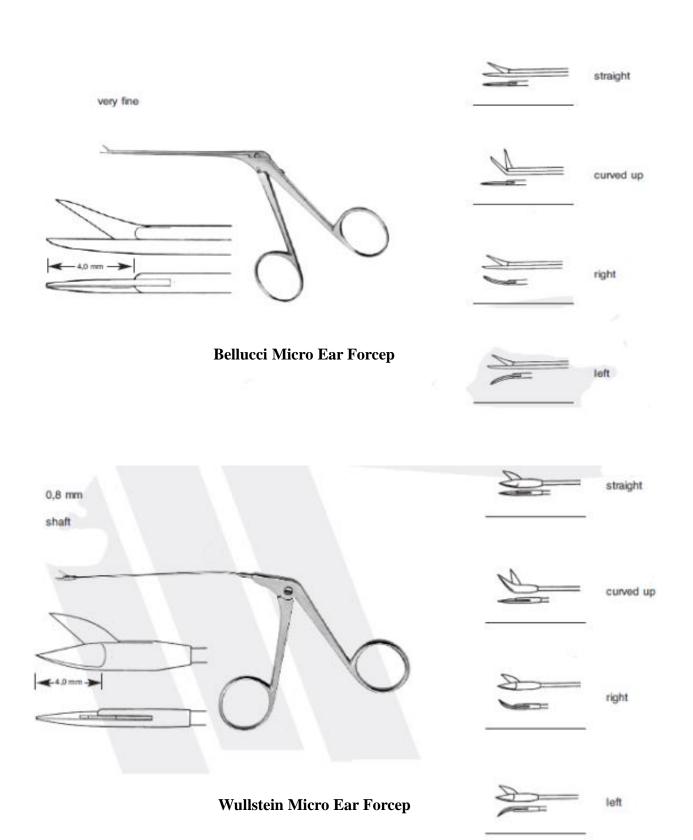

SKIN HOOKLETS, ELEVATORS

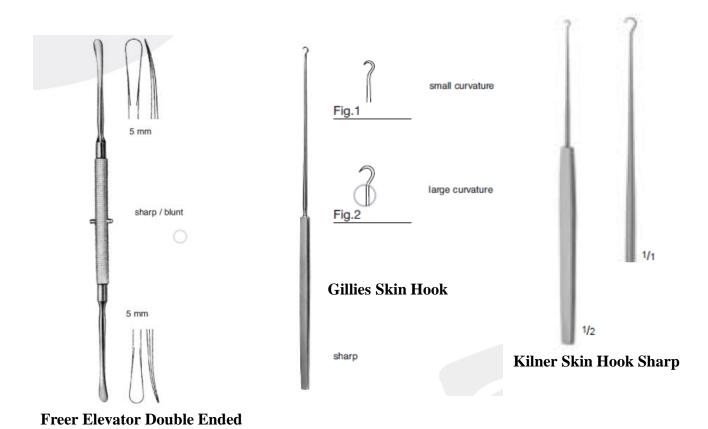

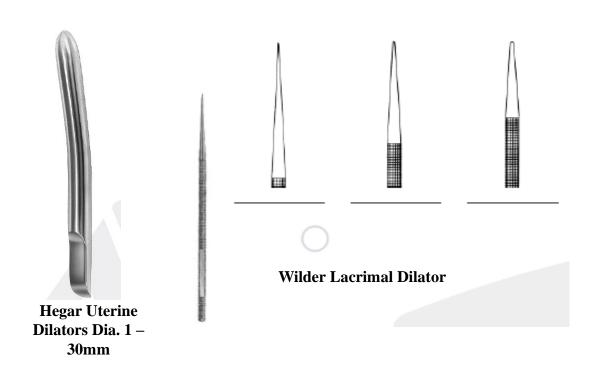

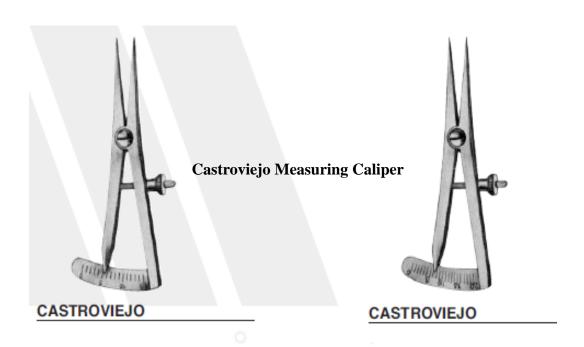

Measuring Range: 15 mm Measuring Range: 20 mm



Iron Intern® Single Arm Retractor

SET-UP

ARH017, ARH022, ARHD22, ARHJ22, ARH026, ARHJ26



Iron Intern® Single Arm #ARH022-ARM

3/D Swinger Rail Clamp #3DSREG







Iron Intern® comes in single and double arms.

Automatic Retractor Holder - ARH 22°

Standard Iron Intern® Single pole with

3/D swinger (U.S Patents). Vertical pole is 22" long. Iron Intern® is a steady retractor always applying the right pressure, enabling you to position the retractor to your exact requirements.

Automatic Retractor Holder - ARH 17"

Iron Intern® Single pole with 3/D swinger (U.S Patents). Single pole with 3/D swinger (U.S Patents). The vertical pole is 17".

Automatic Retractor Holder - ARH 22" Jumbo

Iron Intern® Single pole with 3/D swinger (U.S Patents). Vertical pole is 22" long.

Jumbo Iron Intern[®] is specially designed for Bariatric procedures. Shoulder and elbow joints are larger. Larger joints make the instruments stronger and more resistant to applied pressure.

Iron Intern[®] lends an extra hand in surgery.

Hydra Double Arm

Hydra 22" Iron Intern® (U.S. Patent) double arm allows the surgeon to use two arm holders instead of one. Includes:

- · Hydra Pole
- · 3/D Swinger
- Hydra Yoke

Hydra Jumbo Double Arm


Jumbo Hydra 22" Iron Intern® is especially designed for Bariatric procedures. Larger joints make the instrument stronger and more resistant for applied pressure.

Hydra Yoke

Hydra Yoke is used in Stieber Rib Grip Kit. Attaching the Hydra Yoke to the Lees Center Platform which is located on the Horizontal Bar, gives a surgeon the ability to use two additional arms with that system.


Cowan Gripper

May be used in many procedures. It is an articulated arm which attaches with a C clamp. It may be attached to not only the Iron Intern® system but to other systems, for example a flat ring surgical system.

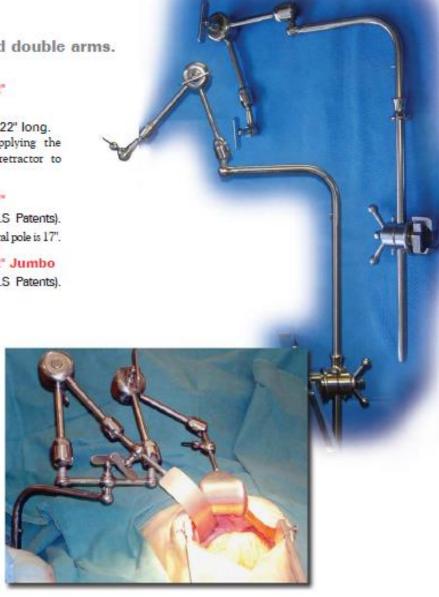



IRON INTERN® FOR LAPAROSCOPIC SURGERY

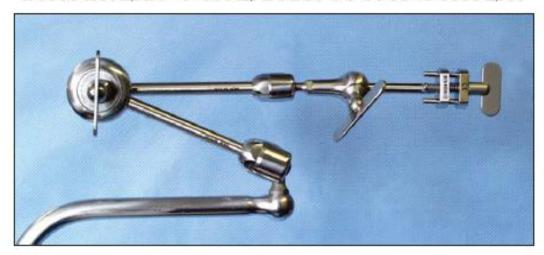

ARH 22" with LC1 Greenstein Gripper

The Iron Intern® is helpful in general, bariatric, and transplantation surgery. The Iron Intern® is adapted for laparoscopic surgery with the use of the Greenstein Grippers and the Nathanson Hook Liver Retractors and other laparoscopic instruments.






ARH 22" with the Nathanson Liver Retractor, medium size

Fields of Surgical Application:

- General
- Laparoscopic
- Bariatric
- Plastic
- Colon-Rectal
- Vascular
- Gynecology
- Transplantation
- Urological
- Pediatric
- Orthopedics
- Thoracic
- Otolaryngology
- · Oral
- Maxillo-Facial
- Veterinary

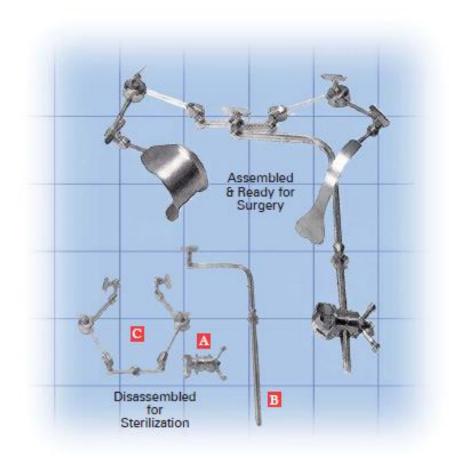

A 3/D Swinger (rail clamp) permits complete versatility of positioning at the site of fixation to the operating table. The 3/D swinger allows the bar to be tilted from side to side 180 degrees. With the Iron Intern*, a simple adjustment of the 3/D Swinger permits rotation of the long bar up to 360 degrees.

B The Vertical Pole allows you to adjust the length of the main support of the articulating mechanical arm by locking it at your required position and tightening the 3/D Swinger on the rail.

Iron Intern® arm functions like a human arm with a shoulder, elbow, and wrist joint. Releasing and locking the wingnut at the elbow joint locks all three. This allows surgeons to adjust the arm to any required position.

The wrist joint holds our large selection of retractors for specific surgery needs. The Universal "Gemma" Connection Clamp (G9) makes using other brand name retractors possible.

